Type 1 Required MODIFICATION INSTRUCTIONS for Service Codes: 3434 and 3425

Kodak X-Omat 3000 RA PROCESSOR Mod No. 1

Kodak X-Omat 5000 RA PROCESSOR Mod No. 3

To help prevent the possibility of film jams in the FIXER/WASH CROSSOVER, due to swelling of the ROLLER, and to install the latest version of software.



Use qualified service personnel to install this modification.

Service Effects: None.
Special Requirements: None.

Serial Numbers: 3000 RA: 0 - 3580

5000 RA: 0 - 756

Special Tools: PORTABLE COMPUTER

INTERFACE CABLE TL-4391

Installation Time: Approximately 0.5 hours

Parts Status: Available from Service Parts Management

Parts Requirements: See the Parts List on Page 2.



HEALTH IMAGING

PLEASE NOTE

The information contained herein is based on the experience and knowledge relating to the subject matter gained by Eastman Kodak Company prior to publication.

No patent license is granted by this information.

Eastman Kodak Company reserves the right to change this information without notice, and makes no warranty, express or implied, with respect to this information. Kodak shall not be liable for any loss or damage, including consequential or special damages, resulting from any use of this information, even if loss or damage is caused by Kodak's negligence or other fault.



This equipment includes parts and assemblies sensitive to damage from electrostatic discharge. Use caution to prevent damage during all service procedures.

Table 1 Mod Kit, No. 1, 9B8376 for the 3000 PROCESSOR includes:

Part No.	Description	Qty.
9B8662	662 "A" ROLLER - FIXER/WASH CROSSOVER	
9B8449	Software Diskette, latest version	
8B7489	Mod Instructions	1
1C0583	Mod Reimbursement Card	

Table 2 Mod Kit, No. 3, 9B8377 for the 5000 PROCESSOR includes:

Part No.	Description	Qty.
9B8662	9B8662 "A" ROLLER - FIXER/WASH CROSSOVER	
9B8247	Software Diskette, latest version	
8B7489	Mod Instructions	1
1C0583	Mod Reimbursement Card	1

Installing the New Software

Note

For more details on loading the new software, see the procedure for Downloading the Operating Software in the Adjustments and Replacements service publication.


- 1. Energize the PROCESSOR, and connect the PORTABLE COMPUTER.
- 2. Insert the DISKETTE from the Mod Kit into the A DRIVE of the PORTABLE COMPUTER.
- 3. At the C:> prompt, type A: and press [Enter].
- 4. At the A:> prompt, type **DOWNLOAD** and press [Enter].
- 5. While the new software is loading, install the new "A" ROLLER. See the procedure below.

Installing the New "A" ROLLER

- 1. Remove the existing "A" ROLLER from the FIXER/WASH CROSSOVER.
 - See Figure 1 for the 3000 RA PROCESSOR on Page 3.
 - See Figure 2 for the 5000 RA PROCESSOR on Page 4.
- 2. Install the new "A" ROLLER from the Mod Kit.

2 01OCT96 - 9B8818

Figure 1 3000 RA PROCESSOR, Replacement of the "A" ROLLER — See Page 4 for the 5000 PROCESSOR

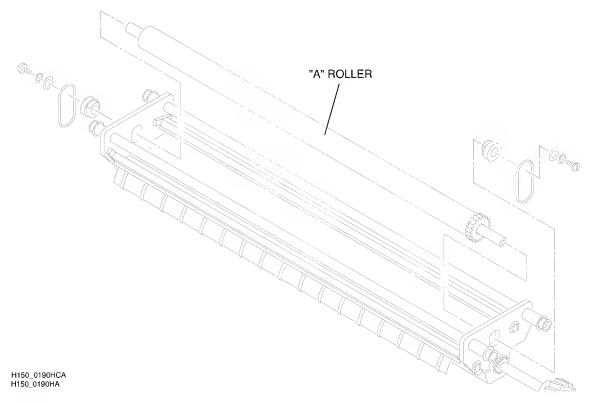


Note

ROLLER 9B8662 does <u>not</u> substitute for 467017 in all positions in the 3000 and 5000 PROCESSORS, only in the "A" position of the FIXER/WASH CROSSOVER.

3. Change the Part No. for the new "A" ROLLER in the Figures in the Parts List, and add the new Part No. to the Numerical and Alphabetical Indexes:

Old	New		3000 RA Processor	5000 RA Processor	
Part No. Part No.		Description	Affected Pages	Affected Pages	
467017	9B8662	"A" ROLLER - FIXER/WASH CROSSOVER	2-5, 3-23, 9-4, and	2-3, 3-23, 9-4, and	
		<u>only</u>	10-7	10-6	

- 4. Circle the "1" on the MODIFICATION LABEL of the $\underline{3000\ PROCESSOR}$ or the "3" on the 5000 PROCESSOR.
- 5. Check for the correct operation of the PROCESSOR.

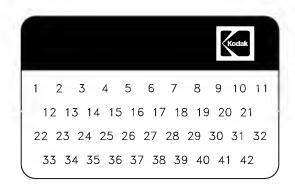


On the MODIFICATION LABEL:

Circle the "1" (3000 PROCESSOR) or the "3" (5000 PROCESSOR)

9B8818 - 01OCT96 3

Figure 2 5000 RA PROCESSOR, Replacement of the "A" ROLLER

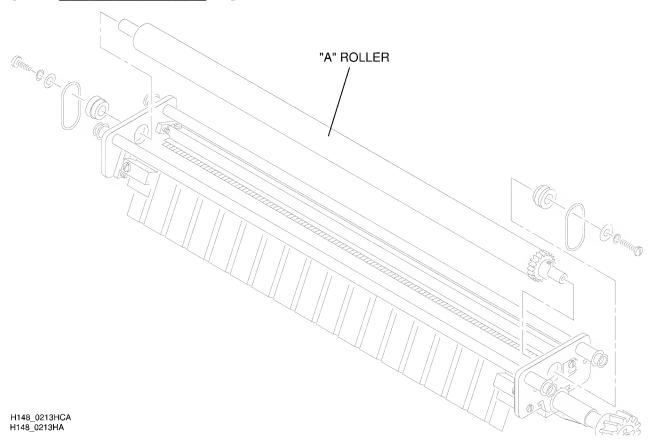

Print Date	Pub. No.	Affected Pages	File Name	Notes
01OCT96	9B8818	All	ma3425_3_01oct96.doc	New Modification Instructions.

Kodak and X-Omat are trademarks.

HEALTH IMAGING EASTMAN KODAK COMPANY • ROCHESTER, N.Y. 14650